Instructions for use

FeelSafe Pro canopy bed BT004

REF BT004

The FeelSafe Pro comes standard with a bed base and an optional mattress topper. The image above is for illustrative purposes only. The canopy, bed base, and mattresses are available separately.

PLEASE NOTE! In urgent situations, the complete canopy bed can be rented immediately. Delivery is 24/7, including Saturdays, Sundays and public holidays. In these situations, please call +31 (0)13 536 38 07.

Please read these instructions carefully to avoid any risks or unwanted bed exit.

KEEP THESE INSTRUCTIONS FOR USE AS A REFERENCE


Human Protection B.V.
Droogdokkeneiland 4
5026 SR Tilburg, The Netherlands
tel. +31 (0)13 536 38 07

Table of contents

Instructions for use Dutch	3
Instructions for use German	37
Instructions for use French	72
Instructions for use English	106

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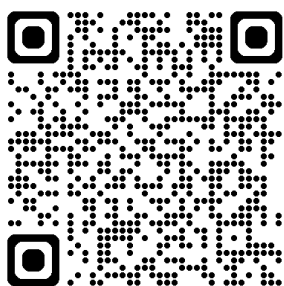FeelSafe Pro tent bed BT004

REF BT004

The FeelSafe Pro is supplied by default on a bed base with an optional top mattress. The above visual is for illustration purposes. The bed tent, bed base, and mattresses are available separately.

ATTENTION! The complete tent bed can be rented last-minute for emergencies.
Delivery 24/7, including on Saturday, Sunday and public holidays.
Call +31 (0)13 536 38 07 in these situations.

Read these instructions carefully to avoid risks or prevent a patient from accidentally exiting the bed.

SAVE THIS MANUAL FOR REFERENCE

Human Protection B.V.
Droogdokkeneiland 4
5026 SR Tilburg, The Netherlands
Phone +31 (0)13 536 38 07

CONTENTS USER MANUAL ENGLISH

Definitions	109
1. Purpose of the device	110
1.1 For whom?	110
1.2 For which purpose?	110
1.3 Indications	110
1.4 Contraindications	110
1.5 Description of the device	111
1.6 Benefits of the FeelSafe Pro tent bed	111
1.7 Classification	111
1.8 Traceability	111
1.9 Reporting defects, complaints or incidents	111
2. Procedure for using the FeelSafe Pro tent bed	112
2.1 Before using the FeelSafe Pro tent bed	112
2.2 First use	112
2.3 During use of the FeelSafe Pro tent bed	112
2.4 End of use	112
2.5 Policy for daily inspection and the patient's freedom	112
2.6 Adjustments to the FeelSafe Pro tent bed	113
2.7 Checklist for first use	113
2.8 Training for care providers	113
2.9 Location of the bed	114
3. 10 important things to check before leaving a patient unattended	115
4. List of parts of the FeelSafe Pro tent bed	116
5. Points to consider when using the FeelSafe Pro tent bed	117
5.1 Other warnings	117
6. Use of the FeelSafe Pro tent bed	118
6.1 Access panels	118
6.2 Self open set	118
6.3 Quick-release fasteners	119
6.4 Zips and pullers	119
6.5 Drainage openings	120
6.6 Transport handles	120
6.7 Universal fastening clamp	121
6.8 Lockable swivel castors	121
6.9 Comfort mattress and mattress compartment	121
6.10 Top mattress	122
6.11 Opening the tent in case of emergencies	122
7. Transfers	123
7.1 Transfer to the FeelSafe Pro tent bed	123
7.2 Access by care providers	124
7.3 Transfer from the FeelSafe Pro tent bed	124
7.4 Assisting patients in emergencies	124
7.5 Which type of hoist to use in combination with the FeelSafe Pro tent bed	124
8. Cleaning instructions for the FeelSafe Pro tent bed	125
8.1 Cleaning	125
8.2 Disinfection	125
8.3 Preparing to wash the tent	125
8.4 Washing the tent	126
8.5 Attaching the tent after washing	126

9.	Maintenance and repairs	127
9.1	Service and maintenance contracts	127
9.2	Technical service	127
9.3	Monthly visual inspection of the tent	128
9.4	Storage and use	129
9.5	Expected service life of the FeelSafe Pro tent bed	129
9.6	Instructions for safely removing the device	129
9.7	Removal	129
9.8	Use of accessories in combination with the FeelSafe Pro tent bed	130
9.9	FAQ	131
9.10	Technical specifications/dimensions	132
9.11	Replacement parts	133
9.12	Warranty	133
10.	Assembly instructions	134

Definitions

A situation that requires your attention

WARNING! A situation that requires your attention A warning alerts you to a hazard that might lead to injury to yourself or to the patient or damage to the FeelSafe Pro tent bed if not prevented.

The following words in this manual have special meanings:

“BIOHAZARD” “ ” refers to contamination with human blood or tissue, certain bodily fluids, recombinant DNA or pathogens. Consult the patient's doctor if you have any questions about this hazard.

“CARE PLAN” refers to the patient's personal care plan, as prepared by the medical team (see below).

“ENTRAPMENT” indicates a situation in which the patient is entrapped, confined, or entangled in a confined space. This may give rise to suffocation, compression, or lack of proper blood flow, resulting in serious injuries.

“MEDICAL TEAM” refers to the interdisciplinary care team that draws up the care plan and oversees the patient's care. This team may include doctors, nurses, occupational therapists, physiotherapists and social workers, among others.

Warning

Medical device

Catalogue no.

Please refer to the manual for instructions for use

Manufacturer

Country of manufacture

Serial number

Keep dry

Keep away from sunlight

1. Purpose of the device

The FeelSafe Pro tent bed consists of a tent and a tent frame, and is mounted on a bed base by default. The FeelSafe Pro tent bed has been developed to offer a safe and controlled environment to patients who are at great risk of injury if they fall out of bed or get out of bed without assistance. The FeelSafe Pro tent bed is a passive alternative to direct, physical fixation. The bed is intended for adult patients with a length of at least 146 cm, a weight of 40 to 165 kg, and a Body Mass Index (BMI) equal or greater than 17.

The bed must always be used in the level, flat position for paediatric patients. The bed restricts the patient's freedom. As such, it must be prescribed by a qualified physician or by the authority that is legally competent for prescribing such medical devices. The bed may only be used under direct supervision and direct authority of the medical team.

1.1 For whom?

This medical device is intended for use by medical professionals.

1.2 For which purpose?

The FeelSafe Pro tent bed gives patients a sense of freedom, comfort and security, by offering them freedom of movement within a contained space with fewer external stimuli.

1.3 Indications

The FeelSafe Pro tent bed is a more humane alternative to physical fixation or sedative medication. Not only does the tent bed increase patient safety, but it also promotes a sense of security and protection. Indications for use of the FeelSafe Pro tent bed include:

- Severe osteoporosis;
- A high risk of fractures due to clinical conditions;
- Insufficient muscle and fat tissue and subcutaneous tissue to absorb the impact of a fall;
- Musculoskeletal disorders that cause an abnormal gait or affect the patient's strength or biomechanics;
- Use of anticoagulants;
- A history of serious trauma as a result of a previous fall.
- Reduced consciousness as a result of medication or neurological disorders that affect alertness and balance and/or as a result of dementia.
- Visual or peripheral nerve disorders that affect consciousness or the body's ability to sense its own position and movements.

Every patient must have an approved CARE PLAN before the FeelSafe Pro tent bed can be used.

The FeelSafe Pro tent bed must be used in accordance with patient care standards and the policy of a legally registered medical facility. All applicable local, regional, national and European directives must also be complied with.

1.4 Contraindications

In the event of one or more of the following contraindications, do not use the FeelSafe Pro tent bed but contact the attending physician for advice.

- Patients who have been diagnosed with a condition that may cause aggressive, violent, or self-destructive behaviour. This may result in self-harm and/or damage to the FeelSafe Pro tent bed. When patients thrash around violently, the FeelSafe Pro tent bed can tip over or move to a position in which the patient may become injured;
- The presence of several lines, drains and/or tubes;
- The need for acute respiratory or cardiac support;
- Cardiogenic or septic shock;
- Claustrophobia;
- PICA eating disorder;
- Suicide risk;
- The patient is smaller than 146 cm, weighs less than 40 kg or more than 165 kg, or has a Body Mass Index (BMI) that is less than 17.
- The use of the tent bed in Category I (critical care) spaces, such as an ICU.
- Use by persons with an increased risk of suffocation is not recommended. This device is not suitable for use if a healthcare professional has established that breathable materials in the bed are necessary to prevent breathing problems.
- Patients who suffer from restlessness, excessive motor activity, or who thrash around frantically in bed, increasing their risk of falling out of bed or of injury.

Due to the nature of the diagnoses, the activity level, or the uncontrolled movements of patients, the FeelSafe Pro tent bed is not suited for all patients. When a patient remains and/or becomes restless in the bed, they should be immediately moved from the bed and another solution be found.

1.5 Description of the device

- Enclosed mattress in a compartment:
Helps prevent the risk of the patient becoming entrapped under the mattress.

- Tent:
The bed is surrounded by a rectangular grey tent with a volume of two cubic metres. The fabric (polyester with a PU + FR coating) is lightweight and skin-friendly. The tent has four access panels, so care providers can get to the patient. Moreover, it has several openings for intravenous lines and/or the call button, among others. The zips are fitted with quick-release fasteners for the patient's safety.

1.6 Benefits of the FeelSafe Pro tent bed

The FeelSafe Pro tent bed was developed to offer the following benefits:

- A controlled environment for patients;
- A soothing effect for some patients;
- Reduced risk of exiting the bed independently;
- Reduced risk of falls of patients;
- Reduced risk of entrapment of patients under the mattress;
- No use of bed rails.
This prevents patients from becoming entrapped between the bed rails or between the mattress and the bed rails;
- Large opening, providing easy access for care providers to administer care;
- Easy access with a hoist thanks to the bed's construction;
- Accessible from 4 sides in an emergency.

1.7 Classification

The FeelSafe Pro tent bed is a Class 1 certified medical device in accordance with Regulation (EU) 2017/745 on medical devices (MDR).

It bears the CE marking in accordance with Annexes I and II of the MDR. This product is subjected to the established internal processes for product verification and validation as laid down in Annex IV of Regulation (EU) 2017/745.

1.8 Traceability

The FeelSafe Pro tent bed's unique serial number provides valuable information with regard to traceability. The unique serial number of the FeelSafe Pro tent bed can be found on the product label on the bed's foot end. Always use this number when submitting reports, complaints or incident reports.

1.9 Defects, complaints or incidents

Always report defects, complaints or incidents directly to us.

Human Protection B.V.
Droogdokkeneiland 4
5026 SR Tilburg, The Netherlands
info@humanprotection.nl
+31 (0)13 536 38 07

2. Procedure for using the FeelSafe Pro tent bed

2.1 Before using the FeelSafe Pro tent bed

- The FeelSafe Pro tent bed must be prescribed by a physician or by the authority that is legally competent for prescribing such medical devices;
- All care providers must receive comprehensive training on how to use the bed prior to using the FeelSafe Pro;
- Make sure that you are fully briefed on the patient's care plan and that you know how often and how rigorously you must check the patient;
- Follow all instructions, warnings, and directions in this manual and all labels on the FeelSafe Pro tent bed.

2.2 First use

- Check the patient several times during the first hours to establish whether the FeelSafe Pro tent bed fulfils the patient's needs. Check whether the patient is comfortable and responds well to the FeelSafe Pro tent bed;
- Human Protection B.V. recommends in-room checks during the first two hours and uninterrupted auditory or visual inspections during the following two hours.

2.3 During use of the FeelSafe Pro tent bed

- Always review the '10-point checklist' on page 115 before leaving the patient unattended;
- Check the patient carefully to establish whether the FeelSafe Pro tent bed fulfils the patient's needs. Check whether the patient is comfortable and responds well to the FeelSafe Pro tent bed;
- Re-assess the necessity for use of the FeelSafe Pro tent bed when the patient's mental or physical condition changes. Only use the FeelSafe Pro tent bed when there is no other less restrictive alternative. Always use the least restrictive device for the shortest possible period of time, taking into account the possibility of accessibility from within (§6.2).
- Continually assess whether the patient responds well to the FeelSafe Pro. If there are signs that the patient does not respond well to the FeelSafe Pro, stop using the device and discuss with the medical team which alternatives are available.

2.4 End of use

- The patient must be assessed and checked before transfer to an ordinary bed or to a hospital bed.
- Check the medical certificate or the care plan or consult the medical team or policy of the institution for instructions about transferring the patient.

WARNING!

NEVER use the FeelSafe Pro tent bed if the tent, access panels or zips are damaged. ALWAYS check the tent and zips before leaving the patient unattended. If you do not heed this warning, the patient can accidentally exit the bed, which may result in a fall and injury.

Importance of patient checks. NEVER leave a patient without proper monitoring. If you fail to check the patient's condition in accordance with the safety policy of your institution, the medical certificate or the CARE PLAN, this can have serious consequences.

Stop using the bed immediately if the patient becomes or remains restless during use of the bed.

2.5 Policy for daily inspection and the patient's freedom

Check the patient's condition in accordance with the safety plan of your institution, the medical certificate and/or the CARE PLAN. Explain to the patient and the patient's family how the call button works.

Follow the instructions in the CARE PLAN with regard to offering the patient daily moments of freedom out of the bed. Take the patient's needs in terms of food, exercise, walks, leisure, toilet visits, and care into account.

2.6 Adjustments to the FeelSafe Pro tent bed

WARNING! NEVER make adjustments to the FeelSafe Pro tent bed for any reason whatsoever.

- ⚠ NEVER add bed rails or boards at the foot or head end of the FeelSafe Pro tent bed. This increases the risk of ENTRAPMENT and serious injury for agitated or aggressive patients who come into contact with these parts;
- ⚠ NEVER use an assistive device, such as a cable tie or padlock to secure the zips on the panels. Always use the quick-release fasteners supplied with the bed. You must ALWAYS be able to get to the patient in case of an emergency;
- ⚠ Do NOT add any accessories to the FeelSafe Pro tent bed other than the ones recommended by Human Protection (see page 130 for a list of accessories).

2.7 Checklist before first use

ALWAYS check whether the FeelSafe Pro tent bed works as it should. Perform a complete visual inspection of the bed before transferring patients to the bed for the first time. Check that:

- The bed is clean and not soiled;
- The bed has been properly assembled and all parts are in good condition;
- The tent frame is not bent or broken;
- All zips can be opened and closed easily and securely and there are no gaps or openings when you apply pressure to the zip;
- The metal tent frame is fully protected and all foam parts are in place and undamaged;
- The tent is not torn and the panels do not show signs of wear, broken seams, or thinning fabric;
- The power cord (if used) is neatly stored so there is no risk of tripping.

NEVER use the FeelSafe Pro tent bed if a check of the above points reveals that the FeelSafe Pro tent bed is not fully functional and/or is damaged. Call us at +31 (0)13 536 38 07.

2.8 Training for care providers

Before use of the FeelSafe Pro tent bed, make sure that all care providers have been fully trained in all applicable policy rules of the hospital with regard to restrictions on freedom and all relevant Dutch and European legislation on involuntary/compulsory treatment. Read and follow the:

- Instructions for use in this manual;
- Warnings and directions in this manual;
- Labels on the FeelSafe Pro tent bed.

2.9 Location of the bed

Make sure that the room is large enough to accommodate the FeelSafe Pro tent bed.

- The space must have a 220V wall socket, within easy reach of the bed.
- In emergencies, you must have enough space to get to the bed easily and remove the patient from the bed;
- Install the FeelSafe Pro tent bed in such a way that patients cannot injure or propel themselves off walls, furniture, or other objects, even when the bed is raised or lowered to a different position.
- Make sure the power cord is not in the way and is not a tripping hazard.
- The space under the bed must remain clear.

Smoking and high temperatures.

- Prohibit smoking in the space where you have installed the FeelSafe Pro tent bed.
- Never expose the FeelSafe Pro tent bed and bed linens to open flames, smoking materials or heat sources. Always keep fan, convector and other mobile heaters away from the bed.

3. 10-point checklist before leaving the patient unattended

Make sure that you are fully briefed on the patient's CARE PLAN or medical certificate and that you know how often and how rigorously you must check the patient.

1. **Tent and mesh** Check that there are no signs of wear and tear, holes or thinning fabric in the nylon panels or mesh. Check that the mesh is visibly taut throughout and does not hang loosely. Check that the tent has been correctly attached to the frame of the bed base and that the metal frame is fully protected by the foam rolls. Check that these are in good condition.
2. **Mattress.** The mattress must always be in the fully enclosed mattress compartment. NEVER place a mattress in the patient space to prevent risk of suffocation. Only the optional Human Protection 5 cm top mattress may be used.
3. **Hygiene.** Make sure the bed is clean and does not contain clutter or objects that may increase the risk of injury or suffocation.
4. **Height and operation of the bed.** The bed must always be in the lowest position, except when care is administered or during transfers. Never leave the bed raised with an unattended resident in it. That way you reduce the risk that the bed tips over when the patient becomes agitated or aggressive. NEVER leave the remote control of the bed base in the patient space when the patient is left unattended.
5. **Lockable swivel castors.** To prevent the risk of unwanted movement of the bed, all four lockable castors on the bed frame must be fully locked when the patient is in the bed.
6. **Zips and quick-release fasteners.** Test whether all zips can be easily and securely opened and closed and there are no gaps or openings when pressure is applied to the zip. Test all zips along their entire length by pressing on the panel or on the tent next to the zip. That way you can check that they are properly closed. ALWAYS use the zip pullers and attach the quick-release fasteners to all the access panels, depending on the degree of freedom you want to provide with the accessibility option (see §6.2)
7. **Location.** Install the FeelSafe Pro tent bed in such a way that patients cannot injure or propel themselves off walls, furniture, or other objects when in the tent.
8. **Drainage openings.** Ensure that all zips are fully closed when the openings are not in use. When using the openings, check that the flow is not obstructed.
9. **Call button.** If the CARE PLAN provides for use of a call button by the patient, you must ensure that the patient can reach this call button easily.
10. **Bed surface position.** Leave the lying surface in a completely flat/horizontal position to avoid creating a gap that the patient can crawl into or behind or in which the patient can become entrapped.

4. List of parts of the FeelSafe Pro tent bed

1. Mattress compartment. Space for a 12-centimetre thick mattress.
2. Tent. The nylon fabric is lightweight, skin-friendly, and water-repellent.
3. Two side panels. One panel on each side with one zip for easy access to the patient.
4. Attachment ring to attach the manual to the bed.
5. Access panels at the head end (5a) and foot end (5b). Provide easy access to the bed's head and foot end. The panels have one zip.
6. Quick-release fasteners. Increase patient safety because the quick-release fasteners can be attached to the zip pullers.
7. Foam protective material. The frame to which the tent has been attached is fitted with 10 protective foam parts to reduce the risk of injury after contact with the frame.
8. Velcro. Attach the side panels to the tent's top panel when the side panels are opened.
9. Drainage openings. The tent has an opening on both sides of the foot end through which IVs and lines can be inserted.
10. Storage pouch. The storage pouch contains the instructions for use and the self-open set (see §6.2).
11. Inner pull cord. For additional freedom, attach the inner pull cord to the zip puller on the inside of the FeelSafe Pro. That way, the tent can be opened from within. This pull cord can be removed. This means the FeelSafe Pro can no longer be opened from within.

5. Points to consider when using the FeelSafe Pro tent bed

The FeelSafe Pro tent bed restricts the patient's freedom. The same rules apply as for other measures that restrict freedom. Also refer to your institution's protocol on these measures.

The use of freedom-restricting measures is never 100% safe; continuous monitoring is required. Patients' behaviour can change from one moment to the next.

Do not leave any clutter or objects in the bed that could increase the risk of self-harm or suffocation.

If a patient is at risk of suffocation, keep the mattress flat with the head end down. If the patient's care plan allows for a call button, check that the patient has access to it.

When leaving the room, ALWAYS check that:

- The bed frame is in the lowest position;
- The brake of the bed frame is engaged;
- All zips and buckles are closed, unless the care plan and/or accessibility feature (§6.2) explicitly allows for an alternative;
- The bed control panel is switched off.
- The bed must always be in the lowest position, except when care is administered or during transfers

5.1 Other warnings

WARNING!

Always consult the instruction for use of the bed base to check how to safely use the bed base.

Do not perform any medical procedures on the bed.

6. Use of the FeelSafe Pro tent bed

6.1 Access panels

The FeelSafe Pro tent bed has four access panels with zips: two side panels, one panel at the head end, and one panel at the foot end. Each zip has a puller that you can use to open and close the zip and can be closed with a quick-release fastener.

The side panels can be opened completely, giving you easy access to the patient. There are Velcro strips at the top, which you can use to keep the side panels in place when they are open.

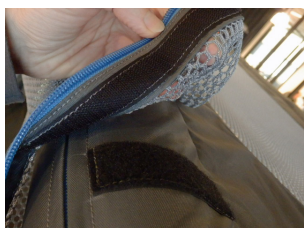
Image: Velcro top of bed

The zip on the right-hand panel has a double slider. The double slider can be used in combination with the self open set to open the FeelSafe Pro from within. See §6.2 for detailed instructions.

Leaving an access panel open can pose a risk, so you must always consider whether this is possible and sensible for the care recipient.

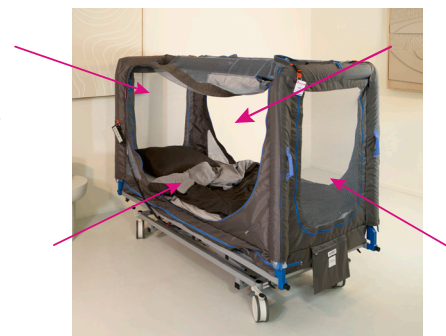
Image: access panels of the FeelSafe Pro tent bed

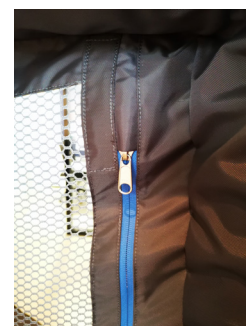
Image: double slider

6.2 Self open set

For additional freedom, you can attach a pull cord to the zip puller on the inside of the FeelSafe Pro. This means the FeelSafe Pro can be opened from within. This pull cord can be removed. This means the FeelSafe Pro can no longer be opened from within.

Attention! When the snap buckles are locked, these will always overrule the option of opening the FeelSafe Pro from within.

6.2.1 Instructions for installation

Step 1 Remove the pull cord from the storage pouch, which also contains the manual.	
Step 2 The pull cord must be attached to the double runner on the inside of the right-hand panel. The pull cord can be installed by looping the black cord through the opening of the puller and then looping the end of the orange section through the opening of the black cord and tightening it.	

6.2.2 Instructions for dismantling the tent-bed

Step 1 The inner cord can be removed by pulling the orange section through the opening of the black cord and then pulling the black cord out of the opening of the zip puller.	
Step 2 Put the inner cord back in the storage pouch.	

6.3 Quick-release fasteners

Each zip has a quick-release fastener. Push the two sides together and pull the clip to open the quick-release fasteners. Push the two sides together to close the quick-release fasteners. You will hear a click when the fasteners are closed. The quick-release fasteners are an additional security to prevent the patient from opening the tent bed themselves.

ALWAYS pull the clip after closing the quick-release fastener to check that the quick-release fastener is securely closed (see visual).

Quick-release fasteners that are damaged or won't close must be replaced immediately!
Call Human Protection B.V. at +31(0)13 536 38 07.


Before leaving the bed unattended, check that the quick-release fasteners on all access panels are securely closed, unless the care plan and/or the use of the accessibility feature (§6.2) allow certain fasteners to remain open. Failure to do so may result in the patient falling or exiting the bed unattended. This could result in serious injury.


6.4 Zips and pullers

- Inspect the zips for wavy or misaligned sections. Test whether all zips can be easily and securely opened and closed and there are no gaps or openings when pressure is applied to the zip.
- If a slider is damaged (see visual below), the zip will no longer close properly.

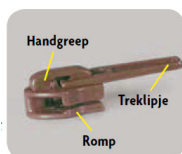Undamaged runner

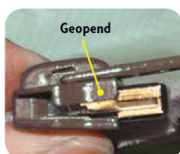Damaged runner

- ALWAYS use the zip pullers to open and close panels.
NEVER pull the panels open: this can damage the fabric and the zips.
- NEVER use the bed if a slider is bent open or if the zip puller is missing or broken.
For repairs, call Human Protection B.V. at +31 (0)13 536 38 07.
- NEVER use the bed when a zip does not work as it should, is damaged, or incorrectly closed.

WARNING! Never leave the patient's bed unattended unless all the zips of the panels have been properly closed, unless the care plan and/or the use of the accessibility feature (§6.2) allow this. ALWAYS test all zips along their entire length by pressing on the panel next to the zip. That way you can check that they are properly closed and whether the access panel can be opened by exerting pressure.

6.5 Drainage openings

Open the zip of an opening and insert an infusion or line through the opening into the patient space. Then close the zip of the opening.

- Check that the line is connected and that the zip is closed as far as possible. Ensure that line flow is not obstructed.

- NEVER insert an infusion or line through the mesh. ALWAYS use the designated openings for lines (see visual above).
- When raising and lowering the bed, check that the lines are not entrapped.

WARNING! Make sure that all zips are fully closed when the openings are not in use. That way, you ensure that there are no openings within the patient's reach that would be big enough for a hand or foot to pass through. NEVER leave an unused opening open if the patient is left unattended.

6.6 Transport handles

When moving the bed to another location, you must always use the transport handles. Before moving the bed, the following safety rules must be observed:

- The bed's surface must be flat and in a horizontal position;
- Unplug the power cord from the wall socket and secure it with the designated Velcro strap on the tent;
- Check the technical condition of the castors;
- Make sure the bed surface is in the lowest position;
- Unlock the central wheel brakes.

After moving:

- Plug the power cord into the wall socket;
- The power cord must be kept away from moving parts of the bed during use and transport.

Examples of uses include: loading and unloading the tent-bed from a transport vehicle, transporting it from the hospital's storage room to the patient's room, and moving the bed from one room to another.

The handles are located on all corners of the bed.

You can adjust the height of the handles by adjusting the height of the bed. Refer to the instructions for use for the bed base to see how to adjust the height.

You can move the bed by pushing or pulling it.

6.7 Universal fastening clamp

A universal fastening clamp can be attached to all four corners of the FeelSafe Pro tent bed. The clamp is attached where the blue frame becomes thicker, as shown in the visual. By aligning the clamp end with the FeelSafe Pro tent bed, the clamp does not pose a risk for people walking past the tent bed. This also prevents damage to the fabric of the tent bed.

Before using such a clamp, always check your institution's protocol.

6.8 Lockable swivel castors

The FeelSafe Pro tent bed is supplied by default on a bed base with four lockable swivel castors. Always lock the swivel castors before transferring patients to reduce the risk of unexpected movement of the bed.

Make sure you lock all swivel castors while administering care to the patient and while transferring patients to and from the bed.

6.9 Comfort mattress and mattress compartment

The FeelSafe Pro tent bed has a closed mattress compartment (see visual). The compartment can accommodate a mattress with these dimensions: 90 cm x 200 cm x 12 cm.

The mattress comes with an incontinence cover with a zip. You should NOT remove this cover. The FeelSafe Pro tent bed is NOT suited for air mattresses.

You can reduce the risk of ENTRAPMENT or suffocation as follows:

- NEVER use a mattress with dimensions other than those for which the compartment was designed. If you use a smaller mattress, gaps will be created between the tent and the mattress. If you use a larger mattress, the fabric will be stretched, damaging the zip.
- NEVER open or close the mattress compartment when the head or knee section of the bed is RAISED. This can damage the fabric and zip.
- NEVER use another mattress than the top mattress supplied by Human Protection B.V. (max. 5 cm thick, 90 x 200 cm) in the patient space.

6.10 Top mattress

The top mattress can be used in the patient space in the mattress compartment to make the bed. Always check the care plan to determine whether the top mattress can be used safely.

The top mattress is never permanently secured in the patient's space, so always check whether the use of the top mattress is appropriate for the patient. Patients who are restless sleepers or who lie awake for long periods of time may end up entrapped under the top mattress. Monitoring is important in order to determine whether the top mattress can be used safely.

Never use another top mattress than the 5-centimetre thick top mattress supplied by Human Protection. The use of a different mattress may increase the risk of suffocation.

6.11 Opening the tent in case of emergencies

In emergencies, always try to open the tent normally first, using the zip on the side panel. If one zip fails to open, you can always revert to one of the other side panels. If these side panels are difficult to access or the zips fail to open, you can use (dressing) scissors to cut open the mesh of the FeelSafe Pro tent bed.

7. Transfers

7.1 Transfer to the FeelSafe Pro tent bed

To ensure a safe transfer and care for patients in the FeelSafe Pro tent bed, a number of factors must be taken into account, including the patient's height and weight, medical condition, limitations, and mobility.

Only use a hoist to transfer patients after checking that there is enough space under the frame of the FeelSafe Pro tent bed to safely use the hoist.

The dimensions of the access panel on the opening side are 110 x 178 cm (H x W) (see visual below). Because the top horizontal bar is positioned towards the centre, the opening provides sufficient space to easily transfer patients to and from the FeelSafe Pro tent bed with most hoists.

Make sure to follow the policy of your institution, the care plan, and the medical certificate to ensure you fulfil the needs of individual patients.

1. Check the area around the bed for tripping hazards. Make sure the transfer can be carried out safely and easily.
2. Check that all swivel castors are locked to reduce the risk of unexpected movements of the bed. Press down the brake to lock the swivel castors.
3. Make sure the bed is flat. The FeelSafe Pro tent bed can be raised up to 95 cm to facilitate patient transfers.
4. Undo the quick-release fasteners and unzip the side panel so the patient can get into the bed. Secure the panel to the top of the bed with the Velcro strips.
5. You can now move the patient to a safe position on the mattress. Depending on the patient's height and weight, you may need additional people for this.
6. If the patient's medical certificate or Care Plan provide for IVs, drains, or monitoring lines or a call button, you must connect these correctly through the designated openings.

WARNING!

If lines or IVs are used for agitated or aggressive patients, then their use in combination with the FeelSafe Pro tent bed must be approved by the doctor or the MEDICAL TEAM.

There is a risk of serious injury if the patient becomes entangled or if the patient pulls out an IV or a monitoring line. Consult with the doctor or MEDICAL TEAM or check the institution's policy for instructions on additional monitoring of such patients.

7. When you are ready to close the tent, lower the side panel from the top of the tent by undoing the Velcro strips. Close the side panel using the zip pullers. Secure the slider with the quick-release fastener. Test all zips along their entire length. Press the panel next to the zip and check that they are securely closed. In addition, check that the access panel cannot be opened by exerting pressure.
8. Read and follow all the steps in the '10-point checklist' on page 115 and ensure that the bed is in the lowest position before leaving the patient unattended.

7.2 Access by care providers

If you want to unzip the tent at the patient's head end, first check that the mattress is FLAT and that the patient is at a safe distance from the panel.

1. Check the area around the bed for tripping hazards.
2. Check that all swivel castors are locked to prevent the bed from moving unexpectedly. Press down the brake to lock the swivel castors.
3. Undo the quick-release fasteners and open the side panel to access the bed. Check that the patient is in a safe position and cannot fall.
4. Secure the panel to the top of the tent using the Velcro fasteners.
5. Check that the patient is not at risk of falling out of the bed when administering care to the patient and when you are ready to open the tent. Carefully move the patient to a safe position on the mattress.
6. Close the access panel on the side and secure it with the quick-release fastener. Check that the bed is in the lowest position.
7. Before leaving the patient unattended, read and follow all the steps in the "10-point checklist" (see page 115).

7.3 Transfer from the FeelSafe Pro tent bed

1. Check the area around the bed for tripping hazards. Make sure that you can perform the transfer safely and easily.
2. Make sure the bed is in a HORIZONTAL/FLAT/LEVEL position. ALWAYS lock the brakes on all swivel castors.
3. Undo the quick-release fasteners and open the side panel to access the bed. Check that the patient is in a safe position and cannot fall.
4. Secure the panel to the top of the tent using the Velcro strips.
5. Help the patient into a sitting position. If necessary, help the patient to move both legs to the edge of the bed and to sit up. Please note: patients who are lying down may feel unsteady or dizzy when getting up. Be prepared and make sure that you are on hand to provide support.
6. Adjust these steps as needed to move patients who need to be lifted out of bed safely. Depending on the patient's size and weight, you may need several people to move the patient safely. Consult your institution's policy for moving patients safely.

7.4 Assisting patients in emergencies

Release the quick-release fasteners and unzip the side panel to help the patient out of bed quickly. If the accessibility feature (§6.2) is activated, the patient may be able to leave the bed independently. If the zip on the side panel fails to open and the zips on the other three panels are difficult to reach or also fail to open, you can use (dressing) scissors to cut open the mesh in an emergency.

A template for a resuscitation plan is available on request from Human Protection B.V.

As resuscitation policies and procedures vary between institutions and hospitals, the resuscitation plan must be tailored accordingly.

7.5 Which type of hoist to use in combination with the FeelSafe Pro tent bed?

Most types of hoists can be used in combination with the FeelSafe Pro tent bed. The leg height of the hoist should be no more than 13 cm. You can also buy or rent a hoist from Human Protection B.V. Call +31 (0)13 536 38 07 for information.

8. Cleaning instructions for the FeelSafe Pro tent bed

The FeelSafe Pro tent bed can be cleaned according to SRI infection prevention guidelines (SRI = Samenwerkingsverband Richtlijnen Infectiepreventie).

For cleaning:
Please refer to the manual of the bed base for information on how to clean the bed base.

8.1 Cleaning

To clean the tent and mattresses, wipe them with a clean cloth with an all-purpose cleaner dissolved in water (follow the dosage instructions on the packaging) or with a damp microfibre cloth. Do not use abrasives on the tent. After cleaning, wipe down with a clean, dry cloth.

8.2 Disinfection

Disinfection is only necessary in case the tent, bed base and mattresses are contaminated with blood or other bodily fluids or in case of final disinfection of a patient with an infectious disease for which disinfection is indicated (such as Norovirus and particularly resistant micro-organisms).

- Blood and other bodily fluids must be thoroughly removed from all surfaces before applying disinfectants.
- Apply a CTGB-approved disinfectant to the surfaces of the bed base, tent and mattress. Follow the instructions on the packaging of the disinfectant. The contact time of the solution promotes efficient disinfection.
- Wipe down the disinfected surfaces with a clean, dry cloth to remove any excess disinfectant.

Remove the mattress and the incontinence cover if necessary. Follow the instructions on page 126 to remove the mattress when it is very soiled. (BIOHAZARD)

Do NOT remove the foam mattress from the cover to clean or disinfect it. Do NOT clean the foam rolls on the frame. This may damage the foam, meaning it can no longer be . The foam rolls are available separately from Human Protection B.V. Call +31 (0)13 536 38 07.

8.3 Preparing to wash the tent

1. Remove all waste from the FeelSafe Pro tent bed;
2. Remove all documents from the storage pouch at the foot end of the bed;
3. Remove the remote control for the bed from the tent;
4. Remove the comfort mattress from the mattress compartment;
5. Remove the top mattress from the tent;
6. Remove the tent;

8.4 Washing the tent

ATTENTION! ALWAYS wash the tent before using it for a new patient. Clean or disinfect the bed base and mattresses. That way you avoid infections and the transfer of micro-organisms.

Make sure all zips and Velcro strips are closed to avoid damaging the tent while washing it.

The tent of the FeelSafe Pro tent bed should be washed at a temperature up to 30°C to prevent shrinkage and wear. Use of a bactericidal detergent is strongly recommended.

Use an extra-large capacity washing machine (at least 18 kg). Wash one tent at a time, using up to half of the load capacity and the high water level setting. Only add detergent to the laundry load after the drum is filled with water. Spin at up to 500 rpm. After the wash cycle, immediately remove the tent from the machine making sure that all water has been removed from the fabric. Air dry the tent. Make sure the tent is completely dry before using it again.

Important warning: If you choose to wash or disinfect the tent yourself, you must always consult with us in advance about the cleaning agent to be used, the washing machine settings and the cleaning and disinfection process. Without explicit approval of the cleaning and disinfection process you intend to use, no guarantees can be given about the washing result, hygiene safety, or service life of the tent bed.

Human Protection B.V. can wash the tent for you, according to a washing process that was specifically designed for this purpose. Contact Human Protection B.V. on +31 (0)13-536 38 07 for further information.

The tent of the FeelSafe Pro tent bed must be checked for tears, holes and excessive wear every time you clean it. Damaged parts must be repaired or replaced. Call Human Protection B.V. at +31 (0)13 536 38 07.

Warnings and precautions regarding washing

- NEVER use bleach, harsh detergents, or solvents. These products can cause the backing of the fabric to come undone, reducing its water-repellency;
- Avoid damaging the fabric and only use detergents intended for washing machines. Make sure all zips are closed and Velcro strips are attached to prevent damage to the tent.
- Too much heat can cause the backing of the fabric to come undone, reducing its water-repellency;
- Check all Velcro strips for dirt on a weekly basis. Brush the Velcro with a brush with strong bristles to remove dirt.
- Please ensure that the Velcro does not come into contact with the mesh during use or cleaning to prevent damage to the mesh.

8.5 Attaching the tent after washing

You will need 2 people attach the tent to the frame.

1. Place the tent on the bed, ensuring that the storage pouch for the instructions is at the foot end of the bed.
2. Unfasten all zips for the upper rails (on the long and short sides).
3. Start on the long side and wrap the section for the top rail around the foam. Pull the tent sections together and close the zip. Repeat this step on the opposite long side.

Ensure that the teeth of the zip are correctly positioned. If the teeth are not correctly positioned, the zip will not close properly. This may lead to the patient leaving the bed unattended and injury.

When zipping the panels closed, make sure to pull the fabric tightly together to prevent damage to the zips or fabric.

4. Unzip all four corner legs.
5. Zip the two short sides closed at the top.
6. Wrap the tent around the foam at all the corner legs. Now pull the tent sections together and close the zip. Make sure the teeth of the zip are well aligned. Use the Velcro to cover the corners. Repeat the step on all four sides until all corner legs are attached to the foam padding of the tent.
7. Open the zip of the mattress compartment, place the comfort mattress (Dimensions 90x200x12cm) in the cover, and close the zip.
8. Place the mattress topper (dimensions 90x200x5cm) in the tent.
9. Fasten the power cord with the Velcro strap at the head end.
10. Attach the remote control to the designated ring on the tent.
11. Place the manual in the storage compartment at the foot end of the bed.

The FeelSafe Pro tent bed is now ready for use.

9. Maintenance and repairs

9.1 Service and maintenance contracts

Human Protection B.V. has various options for maintenance and repairs and after the warranty has expired. Below is a list of options.

- 1) You perform all repairs and maintenance yourself. After the warranty has expired (and only when used correctly), Human Protection B.V. cannot be held liable for incidents that occur and result in injury and/or damage.
- 2) You opt for a maintenance contract for the complete tent bed, with all maintenance performed by Human Protection B.V. Should you have questions about this or if you wish to conclude contract, please call Human Protection B.V. at +31 (0)13 536 38 07.

Inspection of the tent and frame

Your maintenance staff must inspect the FeelSafe Pro tent bed when the tent has been removed for washing. You must also check, at least once a month or before using the tent-bed for a new patient, whether the tent and frame are intact and can be safely used.

Do you have questions about the condition of the FeelSafe Pro tent bed? Contact Human Protection B.V. at +31 (0)13 536 38 07

9.2 Technical service

Tent repairs

Clean the tent and close all the zips to prevent damage during shipping.

Repairs during the warranty period

Send the clean product to Human Protection B.V.*
Don't forget to send your proof of purchase with the purchase date.

Repair outside the warranty period



1. Send the clean product to Human Protection B.V.
2. We will draw up a quote for the repair of each product.
3. After approval of the quote, our workshop will take care of the repair.

Human Protection B.V. has a 24hr repair service, including weekends/public holidays, with the option of a replacement FeelSafe Pro tent bed. Should you have questions about this or if you wish to use this service, please call Human Protection B.V. at +31 (0)13 536 38 07.

9.3 Monthly visual inspection of the tent

Please check the following at least each month or when using the tent-bed for a new patient:

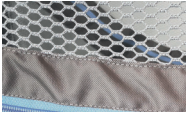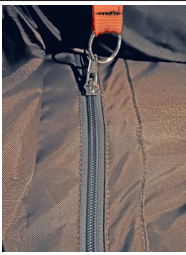		Photo	Clean?		Approved?	
Tent (nylon)	Inside and outside, incl. bed base		Yes	No	Yes	No
Tent: mesh	Check inside and outside of all mesh parts carefully for holes, especially near the edges		Yes	No	Yes	No
Tent: zips and sliders	Check carefully for defects, and the runners in particular. Close all zips, including behind Velcro strips.		Yes	No	Yes	No
Tent: Velcro	There may not be any fluff or threads in the Velcro on top of the tent.		Yes	No	Yes	No
Tent: quick-release fasteners	Check for the presence of defects. Fasten all the buckles.		Yes	No	Yes	No
Top mattress (black)	Check for damage, stains and odours.		Yes	No	Yes	No
Comfort mattress (blue)	Check for damage, stains and odours.		Yes	No	Yes	No
Foam rolls and frame	Present, undamaged and secured.		Yes	No	Yes	No
Check placement of tent on bed base	Check that the tent has been correctly placed on the bed base: head end of the tent aligned with the head end of the bed base. See label on tent and sticker on bed base.		Yes	No	Yes	No
Storage pouch	Check that the following is present: - Instructions for use - Self open set		Yes	No	Yes	No

For inspection of the bed base, please refer to the bed base instructions for use.

9.4 Storage and use

- The FeelSafe Pro tent bed is supplied fully assembled. The technical specifications on page 129 indicate the recommended operating and storage conditions.
- This device was developed for use in normal indoor environments.
- This device can be stored in warehouses at normal ambient temperatures and humidity levels.
- Avoid high humidity levels and excessive moisture, as these can damage the product materials.

9.5 Expected service life of the FeelSafe Pro tent-bed

With average usage, the expected service life of the FeelSafe Pro tent is 5 years or 50 washes, when cleaned and maintained as per these instructions.

ALWAYS follow the cleaning instructions on page 125 and 126 to optimise the bed's service life. You must IMMEDIATELY stop using the bed when a part of it does not work as it should during use. Call Human Protection at +31 (0)13-536 38 07 for a replacement, which can be rented by the day (also call +31 (0)13 536 38 07 in these situations) or to schedule a repair.

WARNING! NEVER use the FeelSafe Pro tent bed when the tent frame is damaged, if the zips cannot be easily or securely opened or closed, if the frame is bent or damaged, if the frame has not been correctly fastened to the bed or if parts are missing. Immediately remove the patient from the bed to avoid serious injury.

9.6 Instructions for safely removing the device

Removal of BIOHAZARDOUS waste. When the fabric or the foam protective material or the mattress of the FeelSafe Pro tent bed becomes heavily contaminated with biohazardous materials, you must remove these parts in accordance with your institution's policy and applicable local, regional, national or European regulations regarding BIOHAZARDOUS materials.

9.7 Removal

Please refer to the instructions for use of the bed base.

9.8 Use of accessories in combination with the FeelSafe Pro tent bed

The following accessories are available for the FeelSafe Pro tent bed:

Item no.	Description	Purpose	Photo
MA-001	Comfort mattress 12 cm	Basis mattress for the compartment of the tent (supplied with the tent bed)	
MA-002	Top mattress 5 cm memory foam	Additional comfort for patients (supplied with the bed)	
VI001	Resuscitation board including cover for attachment to the FeelSafe Pro tent bed	To create a solid base for CPR/ resuscitation	

Only use accessories that are recommended by Human Protection B.V. Other products may pose a risk to the patient or product.

9.9 FAQ

Tent

1. Question: Why is the zip so difficult to open or close?

Answer: There can be several reasons for this.

a. Fabric, bedding or other material is stuck in the zip. You can fix this by pulling back the zip puller gently. Pull the zip in the opposite direction until the zip comes unstuck. Check that the tent is still intact and whether there are holes in the fabric or it has become worn.

b. Check the teeth of the zip to see whether any are broken.

c. Open the zip with the puller. All zips on the access panels feature a plastic label with warnings.

You must replace the label when it is damaged or illegible. Call Human Protection B.V. to request a replacement label.

2. Question: Why does the zip split after zipping?

Answer: The zip slider is damaged. As a result, the zip teeth are not properly aligned when the zip is closed.

WARNING! Do NOT use the FeelSafe Pro tent bed in that case! Call Human Protection B.V. to schedule repairs +31 (0)13 536 38 07 and rent a replacement bed if applicable. You can rent a FeelSafe Pro tent bed by the day. Call +31 (0)13 536 38 07 during and after office hours. We can be reached 24/7, including on Saturday, Sunday and public holidays.

You can avoid this issue by ALWAYS opening the zip with the zip puller. NEVER pull open the access panel. This damages the slider (see visual above). The zip cannot be securely closed as a result.

3. Question: What can I do if the patient scratches the mesh?

Answer: Remove the patient from the tent bed and consult the physician or review the CARE PLAN for information about the right care.

4. Question: How do I clean the tent when it becomes soiled?

Answer: Check pages 125 and 126 for information on how to clean it.

5. Question: How do I remove fluff and/or other material from the Velcro?

Answer: Remove all fluff and materials with a brush with hard bristles. Close all the Velcro before washing the tent to prevent any fluff from ending up in it.

6. Question: Can I use IVs and lines in the FeelSafe Pro tent bed?

Answer: Check "Drainage openings" on page 120.

7. Question: Can I use an inflatable mattress?

Answer: No. The FeelSafe Pro tent bed is intended for use with the supplied mattress. You must always place the mattress in its entirety in the mattress compartment. The zip of the mattress compartment must also always be closed. Check "Comfort Mattress and Mattress Compartment" on page 121.

8. Question: Can I use our own mattresses as a replacement while cleaning the mattress of the FeelSafe Pro tent bed?

Answer: No, the FeelSafe Pro tent bed is intended for use in combination with the supplied mattresses. In exceptional circumstances (e.g., in the event of contamination), another full foam comfort mattress may be temporarily placed in the mattress compartment, but only if its dimensions correspond to those of the mattress compartment (90 x 200 x 12 cm).

WARNING! It is NOT permitted to place any mattress other than the supplied top mattress (5 cm thick) in the tent.

FAQ

1. Question: Where can I find the barcode of the tent of the FeelSafe Pro tent bed?

Answer: You can find the label with the barcode on the pouch with the manual (see visual).

2. Question: Can I use a hoist in combination with the FeelSafe Pro tent bed?

Answer: Yes, check p. 124 for information on the use of hoists in combination with the FeelSafe Pro tent bed. NEVER use a hoist to transfer patients, unless you have checked that there is enough space under the tent frame to safely use the hoist. The dimensions of the access panel on the opening side are 110 x 178 cm (H x W).

For troubleshooting of the bed base, please refer to the bed base instructions for use.

9.10 Technical specifications/dimensions

Dimensions	External (LxW)	approx. 200 x 90 cm
	Lying surface (LxW)	200 x 90 cm
	Clearance under the bed	approx. 13 cm
	Comfort mattress (WxLxH)	90 x 200 x 12 cm
	Tent volume	2 m2
Weight	Tent	8 kg
Conditions for use	Relative humidity	30 - 75%
	Temperature	+10°C to +40°C
Storage conditions	Relative humidity	30 - 75%
	Temperature	+5°C to +45°C

9.11 Replacement parts

Contact Human Protection B.V. to order replacement parts.

Replacement parts for the FeelSafe Pro tent bed

Description	Item no.
Tent, complete	BT015
Cover for top mattress thin	ELM100037C
Cover for comfort mattress	ELM100037A
Top mattress complete 5 cm memory foam	MA-002
Comfort mattress complete	MA-001
Resuscitation board, including cover	VI001
Bolt and rivet nut for the frame	LT-002
Sealed foam padding, per set	OI-FM

In many cases, you can obtain the parts of the bed base directly from us. Please contact us in that case.

9.12 Warranty FeelSafe Pro tent bed

The tent of the FeelSafe Pro tent bed is exclusively designed for use on the specially adapted bed bases PB526II and PB528 and may not be used on any other bed base.

What is covered by the warranty?

Human Protection B.V. Guarantees that the tent for the FeelSafe Pro tent bed does not contain material defects and/or manufacturing faults to the original buyer. If, after investigating the complaint, material defects and/or manufacturing faults are found, we will repair or replace the materials in question free of charge.

How long is the warranty valid?

The warranty is valid for one (1) year from the date of purchase of the tent bed. The invoice date is decisive.

What is NOT covered by the warranty?

Damage to the tent or other parts as a result of:

1. Use in unsuitable conditions or failure to assemble or maintain the product in accordance with our instructions for use and maintenance;
2. Incorrect use (including, but not limited to, damage caused by patients), misuse, modification, exposure to extreme heat or damage caused by natural forces or force majeure (such as flooding or fire);
3. Modification, alteration or repair, unless this work has been carried out by Human Protection B.V.;
4. Consequential or incidental damages, including, but not limited to, loss of income or profits, are not covered.

Maintenance

For maintenance or repair of the tent, always contact Human Protection B.V. Do NOT return the tent to Human Protection B.V. if the tent has not been properly cleaned. Soiled tents cannot be inspected and repaired due to hygiene regulations. In that case, you will be billed the cost for cleaning.

Cost of cleaning and shipping

- We will cover these costs if the product is still covered by the warranty period and if it concerns material defects and/or manufacturing faults.
- After the warranty has expired, these costs will be at your expense.

10. Assembly instructions

Contents Instructions for assembly

Part A - Packing list and parts	40
Part B - Assembly of the frame	41-42
Step 1 - Metal posts for the head end	
Step 2 - Metal posts for the foot end	
Step 3 - Attach the short rails	
Step 4 - Attach the long side rails	
Step 5 - Attach the corner pieces	
Step 6 - Apply foam padding	
Part C - Tent and mattress	43-44
Step 1 - Attach the tent	
Step 2 - Mattress placement	
Step 3 - Secure the power cord	
Step 4 - Remote control of the bed	

Tools for assembly (not included)

1 x spanner, 17 mm

ASSEMBLY INSTRUCTIONS FeelSafe Pro tent bed

We recommend having two people assemble the FeelSafe Pro tent bed to ensure stability during the assembly process and prevent injury.

To simplify assembly, we have broken down the instructions for assembly into three parts:

Part A - Packing list and parts
Part B - Assembly of the frame
Part C - Tent and mattress

* Before assembly, remove all protective packing materials from all parts.
Attention! Do not remove the plastic seal around the foam padding!

Part A - Packing list and parts

Checklist of parts: You must have the following parts. Contact us if any parts are missing. See page 31 for a full list of replacement part numbers.

Fasteners:

- A. Rivet nuts, M12 x 1.5 x 25 (8 pcs)
- B. Bolts, M12 x 1.5 x 25 (8 pcs)
- C. Metal post for right side of head end;
- D. Metal post for left side of head end;
- E. Metal post for right side of foot end;
- F. Metal post for left side of foot end;
- G. Rail for short side of head end (holes to the right of centre);
- H. Rail for short side of foot end (holes to the left of centre);
- I. Long side rails for corner mounting;
- J. Long side rails for centre mounting;
- K. Corner piece (2 pcs);
- L. Foam rubber padding, short (2 pcs);
- M. Foam rubber padding, long (8 pcs).

Part B: Assembly of the frame:

You will need two people to assemble parts 1 through 4.

The frame consists of various metal parts. The U-shaped profiles are for the bed's head and foot end.

Step 1: Metal posts for the head end (see page 42 for a visual of the bed frame)

Assembly step 1:

Insert the metal post for the right side of the head end (part C) on the right side of the head end into the bed frame (see visual). Make sure that the holes at the top of the post are facing inwards.

Assembly step 2:


Install the adjustable clamp lever on the underside of the bed frame. Tighten it securely.

Assembly step 3:

Repeat on the left side of the head end with the metal post for the left side of the head end (part D).

Assembly step 1:

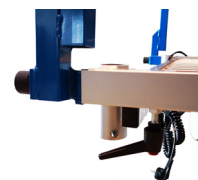Assembly step 2:

Step 2: Metal posts for the foot end

Assembly step 1:

Insert the metal post for the right side of the foot end (part E) into the bed frame on the right side of the foot end. Make sure that the holes at the top of the post are facing inwards.

Assembly step 2:

Install the adjustable clamp lever on the underside of the bed frame.

Assembly step 3:

Repeat on the left side of the foot end with the metal post for the left side of the foot end (part F).

Step 3: Attach the short rails

Assembly step 1:

Ensure that there are two people at the head of the bed and take the rail for the short side of the head end (part G) (holes to the right of centre). Align it between the metal posts on the right and left sides of the head end (parts C and D). Make sure that the holes on the short rails are facing inwards.

Assembly step 2:

Place a rivet nut on an M12 bolt at one end of the perpendicular connection and insert it into the hole. Tighten the bolt with the 17-mm spanner. Repeat this step on both sides of the metal posts.

Assembly step 3:

Repeat steps 1 and 2 at the foot end of the bed and attach the other short rail for the foot end (part H) (holes to the left of centre) to the metal posts on the right and left sides of the foot end of the bed (parts E and F).

Step 4: Attach the long side rails

Assembly step 1:

Repeat this with the long side rails for centre mounting (part I) on the right long side of the bed between the two U-shaped frame parts you have just assembled. Use four M12 bolts with rivet nuts to secure the frame parts. Tighten the bolts with a 17-mm spanner.

Assembly step 2:

Repeat this with the long side rails for centre mounting (part J) in the centre of the bed between the two U-shaped frame parts you have just assembled.

image: step 2

image: fully assembled frame

Step 5: Attach the corner pieces

Assembly step 1:

Place a corner piece (part K) on the two remaining corners. Use four M12 bolts with rivet nuts to secure. Tighten the bolts with a 17-mm spanner.

Step 6: Apply foam padding

Assembly step 1:

Fold open the foam padding (parts L and M) and place them, including the plastic seal, over the frame parts so that the entire frame is covered. Make sure the openings of the tubes face the outside of the frame (see image below).

image: full coverage of the frame

Part C: Tent and mattress

Step 1: Attach the tent

Assembly step 1:

Fold the tent over the bed, ensuring that the storage pouch for the instructions for the FeelSafe Pro tent bed is at the foot end of the bed.

Assembly step 2:

Undo all the zips.

Prevent damage to the zip: pull the halves of the zip towards each other when closing the tent's zips.

Assembly step 3:

Wrap the zip around the foam padding of one of the two long beams at the top, pull the zip halves together and zip it closed. Repeat for the other long beam. Repeat for both short sides. See visuals below.

image: closing the long side

image: closing the short side

Assembly step 4:

Close the zips on all four posts in the same way as in assembly step 3. The zip on the posts extends to the top of the tent on two sides.

Assembly step 5:

Cover the corners with the corner cover which you secure to the Velcro on the posts.

Assembly step 6:

Close all the zips and secure them with the buckle fastener. The buckle fastener must be closed when the bed is left unattended unless the care plan and/or the use of the accessibility feature (§6.2) provide for an alternative.

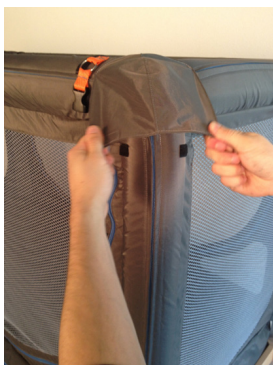
image: covering the corners

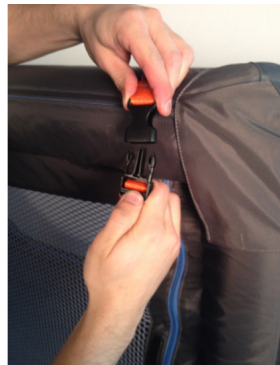
image: closing the buckle fastener

Step 2: Mattress placement

Assembly step 1:
Open the zip of the mattress compartment.

Assembly step 2:
Place the mattress in the mattress compartment.

Only use a mattress with the following dimensions in the mattress compartment: 12 cm x 90 cm x 200 cm. If you use a smaller or larger mattress, the tent of the FeelSafe Pro tent bed may be damaged. The use of a smaller or larger mattress may also result in patient injuries.

Assembly step 3:
Zip up the mattress compartment completely after placing the mattress in it (see visual).

image: zip for closing the mattress compartment

Assembly step 4:
Place the top mattress in the patient space.

Step 5: Secure the power cord
Use the Velco strip at the head end of the tent to secure the power cord during transport.

Step 6: Remote control of the bed
Attach the clip on the back of the remote control to the ring on the head end of the bed.

The FeelSafe Pro tent bed is now fully assembled.
The manual with instructions for use must be stored for reference in the storage space of the FeelSafe Pro tent bed at the foot end.

human
protection